

## Ultrafast digitizers used to identify neutrinos in huge international project

*ADC cards by Spectrum Instrumentation used for largest liquid-neutrino-detector of mankind*

**Grosshansdorf, Germany – 04.09.2024.** Until recently, neutrino particles were thought to have no mass but now the theory is that they actually have mass with a very tiny value as well as coming in three different ‘flavours’ which they can switch back and forth. Often called ghost particles, studying them is very difficult as they typically pass through most normal matter unimpeded and undetected, so specialist detectors have to be built to hunt for them. The latest one is called JUNO, located 750m below ground in Jiangmen, China, and it is made possible by 17 different countries with 730 scientists working in 74 Universities and national laboratories, that have joined forces on this 400-million-euro project. To develop the core part of the detector, the liquid scintillator, ultrafast digitizer cards by Spectrum Instrumentation are used.

JUNO is precisely positioned between eight existing nuclear reactors that provide a source of neutrinos for study. At its heart is a gigantic, highly transparent, acrylic sphere with an inner diameter of 34.5 m, filled with 20,000 tons of a specially developed oil-like substance. This liquid scintillator creates photons when a neutrino interacts with it, and it is encased in a large 35,000 tons water pool. The photons are detected by an array of ca. 45,000 photo multiplier tubes (PMTs) that surround the sphere. The teams at the Technical University of Munich (TUM) and at the Johannes Gutenberg University in Mainz are using M4i.2212 digitizer cards by Spectrum Instrumentation in their high-precision, lab-scale experiments to characterize the liquid scintillators, which require advanced data acquisition. When the JUNO detector is commissioned at the end of 2024, it will be the largest, liquid neutrino detector built by mankind. The detector will drastically enhance our knowledge about interactions and properties of these elusive ghost particles.



*The JUNO main detector 750m underground in a dedicated laboratory. The photo shows the (still empty) water pool of the detector with the central scaffolding. Inside this globe, the acrylic sphere with 34.5m diameter is placed, filled with the liquid scintillator. The white cover just protects sensitive components during mounting.*

### Neutrino detection

The central acrylic sphere contains the liquid scintillator surrounded by a layer of water. Both have to be ultra-pure as the smallest amount of contamination could contain a radioactive material. During construction, workers had to wear two pairs of gloves as even the sweat from a fingerprint could contaminate and ruin the entire project. The detector sits in a specially dug space that is 750 m below ground to shield it from ambient radiation.

When a neutrino interacts with the liquid scintillator (LS), it deposits the energy of the interaction with the molecules of this substance. The enormous light output of the LS (typically > 10.000 Photons / MeV) ensures a precise determination of the deposited energy. It would be highly beneficial if the direction of the incident neutrino could also be reconstructed. Here, the faint but directional Cherenkov light from neutrino’s initial passage through the water is paired to give physicists this information.

---

### Headquarters

Spectrum Instrumentation GmbH, Germany  
Phone: +49 4102-6956-0  
Email: [Info@spec.de](mailto:Info@spec.de)

### US Office

Spectrum Instrumentation Corp., USA  
Phone: (201) 562-1999  
Email: [Sales@spectrum-instrumentation.com](mailto:Sales@spectrum-instrumentation.com)

The aim of the current development of liquid scintillators in Munich and Mainz is to separate the fast but faint Cherenkov light from the dominant scintillation light in order to enable simultaneous energy and directional reconstruction. Therefore, the team under Dr. Hans Steiger constructed several, precision, table-top experiments with enhanced light-collection capabilities and time-resolution.

“We chose the digitizer cards from Spectrum as they provide us with state-of-the-art performance but, unlike rival offerings, they are not expensive or custom creations,” said Dr. Hans Steiger, who is leading the project. “Spectrum’s modular approach means that we could specify exactly what we needed the cards to do so we were not having any compromises or wasting money on unwanted features. I love the fact that they are a standard PCIe product so that we can expand the system on a standard computer chassis as we received more funds. As a university taking part in major long-term international projects, we need to have reliable parts and Spectrum’s five-year warranty gives us piece of mind.”


### JUNO results also push astronomy research

Beyond the work on event reconstruction, the team contribute a calibration project to JUNO. This characterizes the detector material using radioactive gamma and neutron sources where the energy and incident direction are pre-determined. “Our characterizations of the liquid scintillators are only possible because of the ultrafast digitizer cards that enable us to work with timeframes that are measured in picoseconds. Also, the dynamic range of 5V is much better than rivals that are typically 1V which means they can easily manage the 3V pulses in our PMTs that we encounter,” Meishu Lu, a PhD student in the TUM group pointed out. And Manuel Böhles working in Mainz added, “Spectrum has been really supportive in helping us working out the best solutions for our project and helping resolve any issues that we encounter with a phone call directly with one of their engineers. It’s great that they are committed to supporting fundamental research in many universities such as ours.”

The diagram shows the first pulse from the Cherenkov radiation followed by the scintillation signal that gives the energy information. This happens in less than two nanoseconds. By combining this information, the type of particle can be determined and where it came from. This could be from one of the Chinese reactors, the sun, the heart of the earth or deep space. “We have never been able to know exactly where a neutrino was coming from in scintillation detectors before, so this opens up whole new areas of research,” explained Dr. Steiger. “If, for example, a dying star, or so-called supernova, emits large amounts of neutrinos in the sky, we are now able not only to see the neutrinos but also reconstruct with high precision, the point in the sky, where this explosion took place. Effectively we now have a telescope to look into these different neutrino sources to better understand the processes. By detecting light over the entire spectrum, gravitational waves and now also neutrinos with high statistics, energy resolution and directionality start a new era of multi-messenger astronomy.”

### About Spectrum Instrumentation

Spectrum Instrumentation, founded in 1989, uses a unique modular concept to design and produce a wide range of more than 200 digitizers and generator products as PC-cards (PCIe and PXIe) and stand-alone Ethernet units (LXI). In over 30 years, Spectrum has gained customers all around the world, including many A-brand industry-leaders and practically all prestigious universities. The company is headquartered near Hamburg, Germany, known for its 5-year warranty and outstanding support that comes directly from the design engineers. More information about Spectrum can be found at [www.spectrum-instrumentation.com](http://www.spectrum-instrumentation.com)

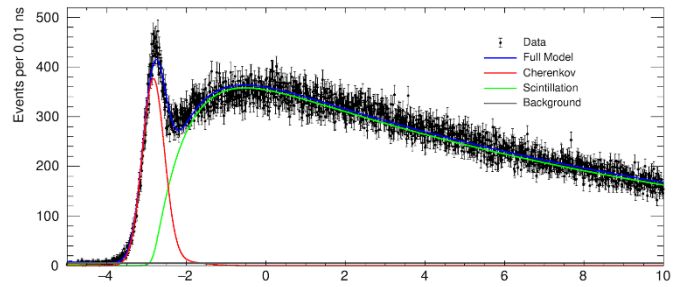

*Typical light emission kinematics for a slow liquid scintillation mixture. The Cherenkov light (red line) in form of a sharp peak in time is followed by the slower scintillation light decay (green line)*



*The M4i.2212-x8 PCIe digitizer by Spectrum Instrumentation with 1.25 GS/s sampling speed on 4 channels.*

### Headquarters

Spectrum Instrumentation GmbH, Germany  
 Phone: +49 4102-6956-0  
 Email: [Info@spec.de](mailto:Info@spec.de)

### US Office

Spectrum Instrumentation Corp., USA  
 Phone: (201) 562-1999  
 Email: [Sales@spectrum-instrumentation.com](mailto:Sales@spectrum-instrumentation.com)